


COMPAT-ONE

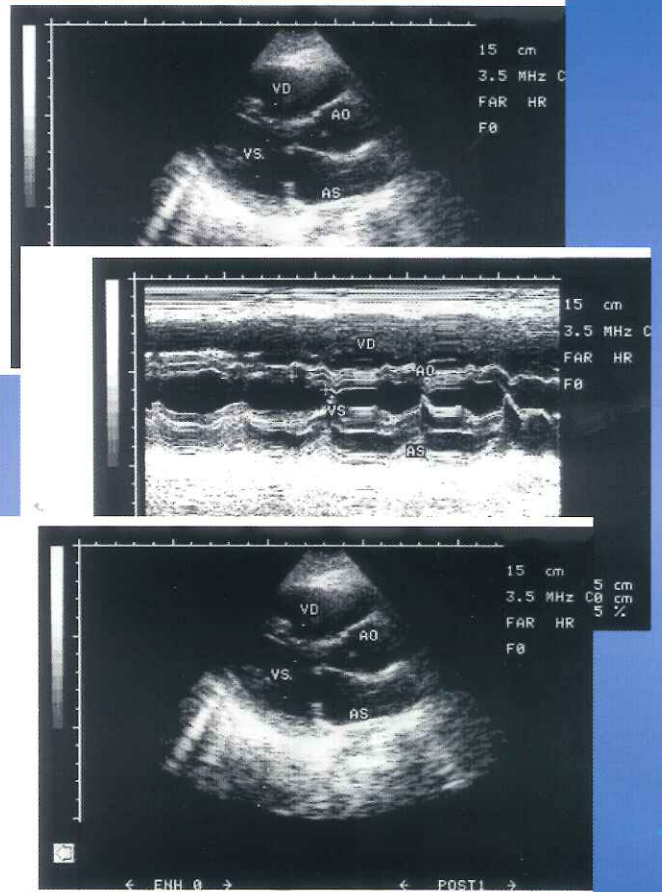

Video printer papers

110S - 110HD - 110HG
KP61B - KP65HM - KP91HG

Compatible with the following printers:

UP-D898MD - UP-X898MD
P93 - P95DE - P95DW

WHY COMPAT-ONE?

- **Excellent performance at unbeatable prices**


This is the paper delivering the best performance with the most profitable conditions; specially developed and manufactured to meet the specific printing requirements for the major brands on the market (Mitsubishi and Sony).

- **Guarantee of durability and stability**

Guarantee of durability and stability under different climatic conditions, like humidity and heat, thanks to the technical features of the chemical composition of Compat-One, particularly the topcoated layer applied on the surface of the paper.

- **High resolution and high image quality**

The structural composition of Compat-One enables, thanks to the function carried out by each single layer, the best performance as far resolution and image quality are concerned.

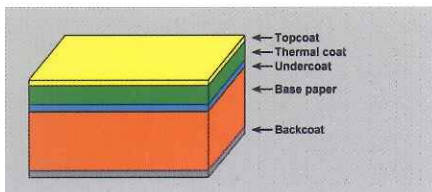

The Topcoat layer

It prevents from deterioration caused by humidity, allows a long conservation in time and extends the service life of the thermal head of the printer.

- **Excellent quality of grayscale reproduction**

The thermal coat, composed by chemical elements specifically designed, reacts when it gets in contact with the thermal printer, changing the color of the paper, ensuring an excellent printing quality and an excellent quality of grayscale reproduction.

The undercoat applied to the base paper guarantees an even and smooth surface, enabling high resolution and high quality, optimizing the sensitivity characteristics of the thermal coat.

- **Better preservation of the printer**

A backcoat layer in the reverse side prevents from accumulation of electrostatic charges, enabling a better preservation of the printer, and from external influences, oils, organic solvents, humidity and water.



Certificate n.
9124.PIR2
EN ISO
13485:2012

CERTIFICATES IN FORCE



Certificate n.
9120.PIRR
ISO
9001:2008



COMPAT-ONE is a brand by Pirrone S.r.l.

Via Balduccio da Pisa 12 - 20139 Milano - Tel. +39.02.52.53.97.55 - Fax. +39.02.57.40.40.77

info@pirrone.com - www.pirrone.com